

Single Use Instruments

Absorbents

Wipes

Woundcare

Patient Identification

First Aid

Cotton Wool

Personal Care


Continence Care

Professional & Personal Healthcare Catalogue





Professional Healthcare

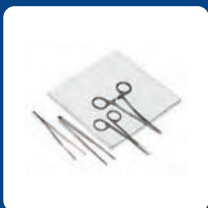

Single Use Instruments



Absorbents



Wipes



Woundcare



Patient Identification

Namelets® ID Bands

Namelets® patient identification bands provide an established method of secure and safe patient identification for any healthcare environment. Easy to use, Namelets® are made from a high quality medical grade polymer that is free from natural rubber latex. This material is also soft and conformable with no sharp edges, minimising the risk of discomfort during wear.

Namelets® are available in two variants. The write-on version has a special scratch resistant panel that can be written on using any normal ballpoint pen. Special pens are not required. The information remains intact through washing or bathing, securely maintaining integrity of the data. Alternatively, an insert version features a card that can either be hand written or computer generated. This is inserted into a pocket that is sealed by the band to prevent water damage or tampering.

The band is fastened using strong plastic stud closures. Namelets® are designed to prevent re-use or transfer between patients. If the band is stretched for removal it becomes visibly deformed, identifying that the patient has tried to remove it. This ensures that patients are identified correctly as it alerts staff to re-confirm all personal details if they see the band is stretched. For safety, bands will also break naturally if sufficient force is applied, to prevent injury to the patient if the band inadvertently becomes entangled. Bands should be cut with scissors when removal is required.

Namelets® are available in adult and child sizes. The child size also fit newborns. The bands are supplied in white or red in write-on style and clear or red in the insert version. Write-on versions are also compatible with most electronic tagging systems.

Red bands are used to alert staff to a particular patient risk, such as an allergy or patient preference. White bands are used as standard for all other patients.

Security

- Strong plastic stud closure
- Cannot be transferred or re-used without being visually detected
- Once written on, information remains intact

Safety

- Band breaks safely in extremis
- Fastening studs are firmly secured and cannot fall out
- Accepts most electronic tagging systems
- Option of red bands, to alert staff to a particular patient risk

Comfort & Ease of Use

- Soft, conformable medical grade polymer for maximum comfort
- No sharp edges to damage delicate skin
- Suitable for any patients
- Available in 2 sizes; child and adult



Product Code	Style	Colour	Size	Unit of Sale
6030	Write-On	White	Adult	100 (1 Box)
6031	Write-On	White	Child	100 (1 Box)
6032	Write-On	Red	Adult	100 (1 Box)
6033	Write-On	Red	Child	100 (1 Box)
6034	Insert	Clear	Adult	100 (1 Box)
6035	Insert	Clear	Child	100 (1 Box)
6040	Insert	White	Child	100 (1 Box)



LATEX FREE



Insert Style

Robinson Healthcare Ltd

Lawn Road
Carlton-in-Lindrick
Worksop
S81 9LB

Tel: +44 (0)1909 735000

Fax: +44 (0)1909 735040

Email: enquiry@robinsonhealthcare.com

Web: www.robinsonhealthcare.com

